



CE 0123



INNOCATH™

IV Catheter with Injection Port

• *Smooth catheterization design*



- Maximizing patient's comfort and caretaker's convenience
- A safe, blood leakage proof IV Catheter
- Ultra sharp triple-facetted bevel needle
- Atraumatic catheter tip design

Smooth Venipuncture with


- Ultra sharp triple-facetted bevel needle
- Atraumatic catheter tip design
- Specially treated catheter tube with excellent gliding characteristic

Needle Free Injection Port

- Enables additional injection without the use of hypodermic needle

Integrated Luer Lock Connector

- For easy and tight closure of the catheter after withdrawal of the needle

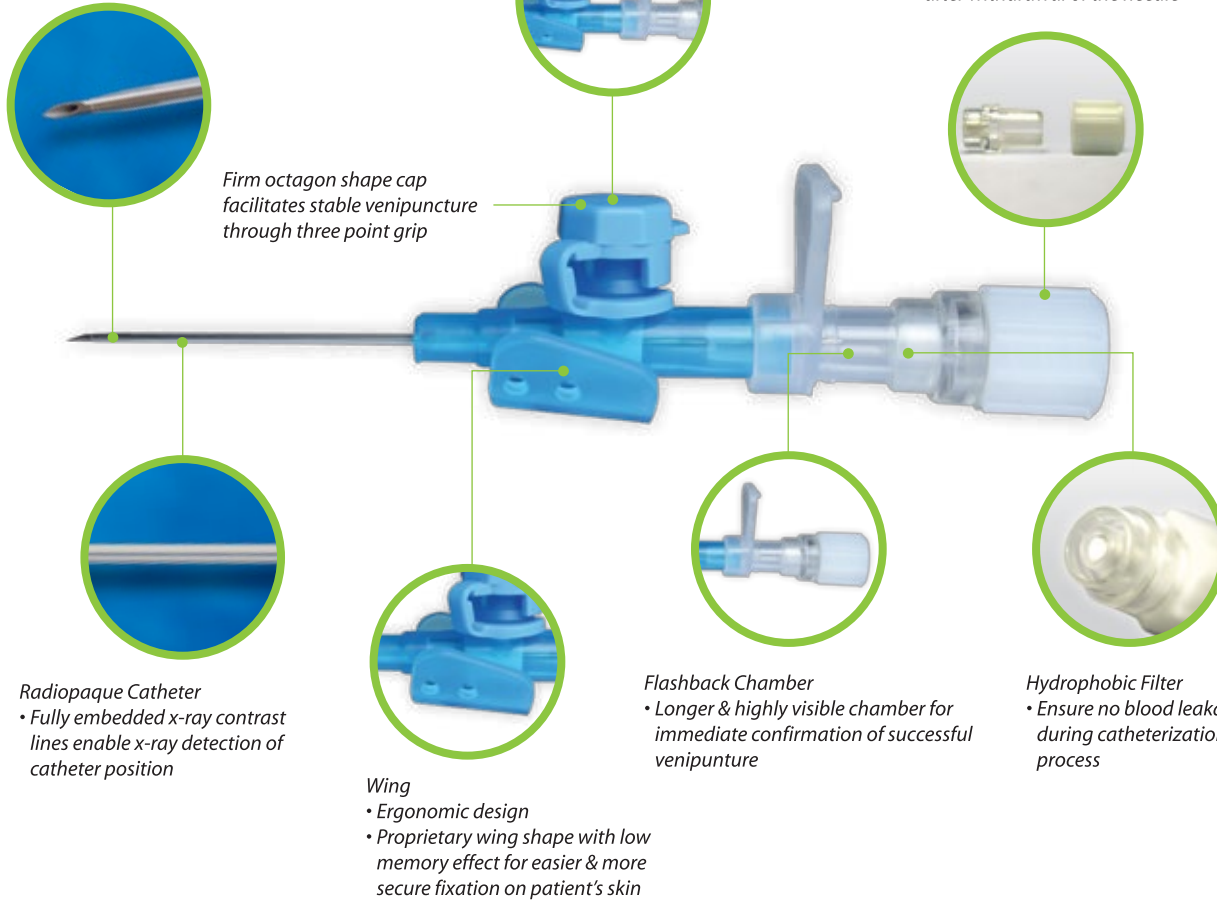

Firm octagon shape cap facilitates stable venipuncture through three point grip

Radiopaque Catheter

- Fully embedded x-ray contrast lines enable x-ray detection of catheter position

Wing

- Ergonomic design
- Proprietary wing shape with low memory effect for easier & more secure fixation on patient's skin

Flashback Chamber

- Longer & highly visible chamber for immediate confirmation of successful venipuncture

Hydrophobic Filter

- Ensure no blood leakage during catheterization process

- Full-range product line•••
- For various clinical indications and therapies
- A wide range of catheter sizes and lengths
- Catheter material made of **FEP**

Size	14G	16G	18G	18G	20G	22G	24G
Colour	Orange	Grey	Green	Green	Pink	Blue	Yellow
Catheter (I.D. mm)	1.7	1.3	0.9	0.9	0.8	0.6	0.5
Catheter (O.D. mm)	2.1	1.7	1.3	1.3	1.1	0.9	0.7
Catheter Length (mm)	45	45	45	32	32	25	19
Water Flow (ml./min.)	300	190	90	98	58	34	21
Qty. per Box (pcs.)	50	50	50	50	50	50	50
Qty. per Carton (pcs.)	500/1000	500/1000	500/1000	500/1000	500/1000	500/1000	500/1000



www.welfordmedical.com

Welford (UK) Ltd
 145-157 St. John Street,
 London, EC1V 4PW,
 United Kingdom.
 T: +44 207 6085 503
 F: +44 203 6425 919
 E: sales@welfordmedical.com

Regional Offices:

Welford Healthcare Sdn Bhd - Malaysia
 T: +603 5567 9336 F: +603 5567 9348

PT Welford Indonesia - Jakarta
 T: +6221 2358 5649 F: +62 21 2358 5601

Welford Vietnam Co. Ltd
 T: +84 4 3 633 1375 F: +84 4 3 633 1401 - Hanoi
 T: +84 8 3 520 2822 F: +84 8 3 823 7840 - HCMC