

SELLA

LINET

Well-Being & Independence



Multifunctional Chair With Excellent Manoeuvrability

Benefits of Sella

For your well-being

- Independence and self-reliance
- High comfort
- Hydraulic lift for better mobilisation



- Reduced physical effort
- Simple examination and easy access



- An attractive design
- Increased efficiency of the hospital ward



Patient


Staff

Management

Excellent manoeuvrability

In combination with the necessary functions, ergonomically suitably selected proportions and a friendly design, the Sella multifunctional chair brings a number of advantages as well as ease for the immobile person and service staff.

Its low weight, great manoeuvrability and sitting comfort make this product an excellent helper in nursing homes, long-term medical institutions, hospitals and also in private care.

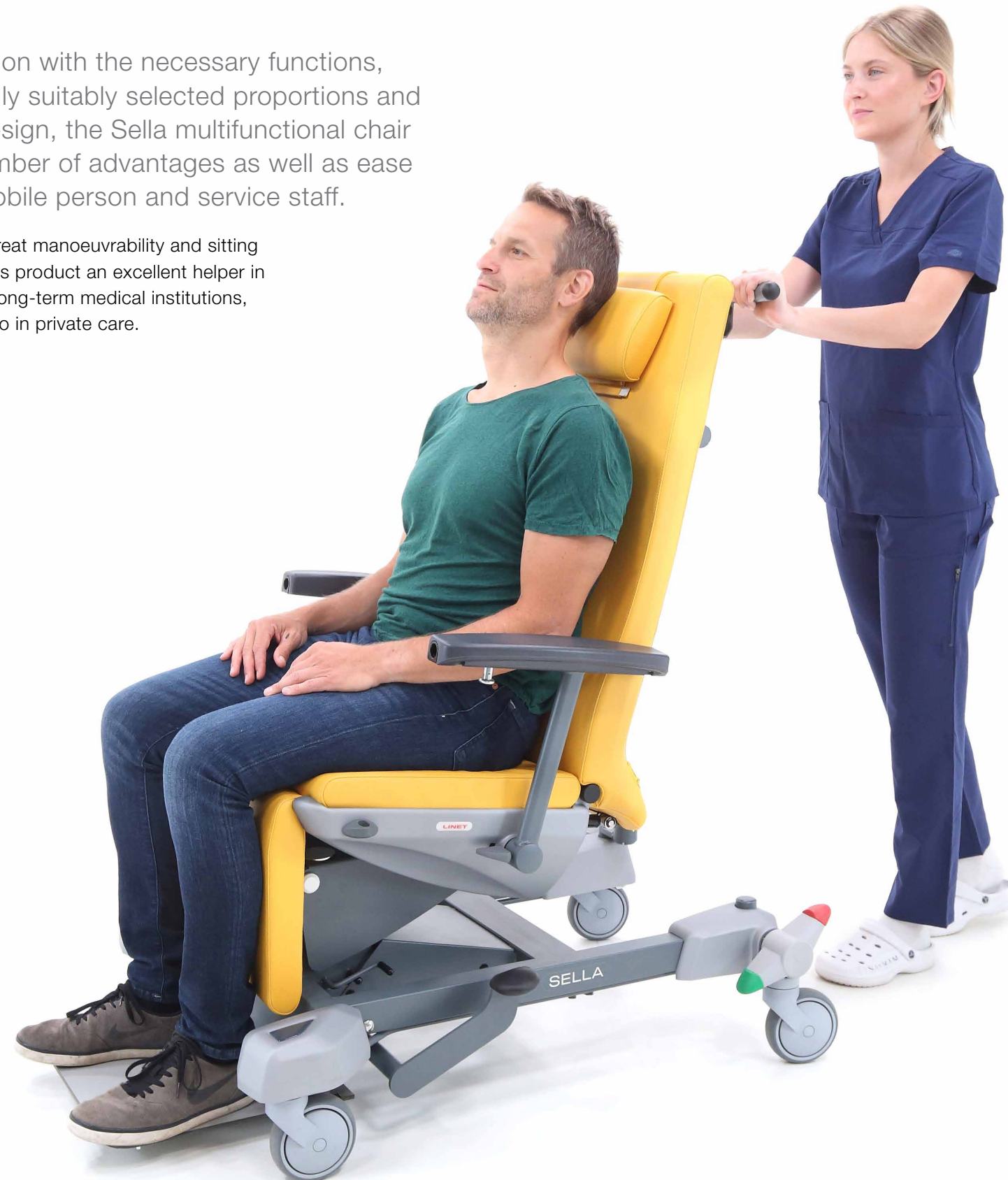

Rotation of the castors

A chair with small castors is more appropriate for comfortable transport in small and narrow spaces. The rotation of the castors around an axis increases the manoeuvrability of the chair in open and small spaces.



Easy transport



Not only does the low weight of the chair contribute to comfortable transport but also to ergonomic control.



The central brake increases patient safety and the simplicity of control for the nurse.


Easy mobilisation

Mobilisation

Hydraulic lift for tipping the seat and adjusting the armrests for arm support completely eliminates the physical exertion of the staff that is necessary for manipulating the patient.



Mobilisation is one of the most physically demanding procedures in healthcare.

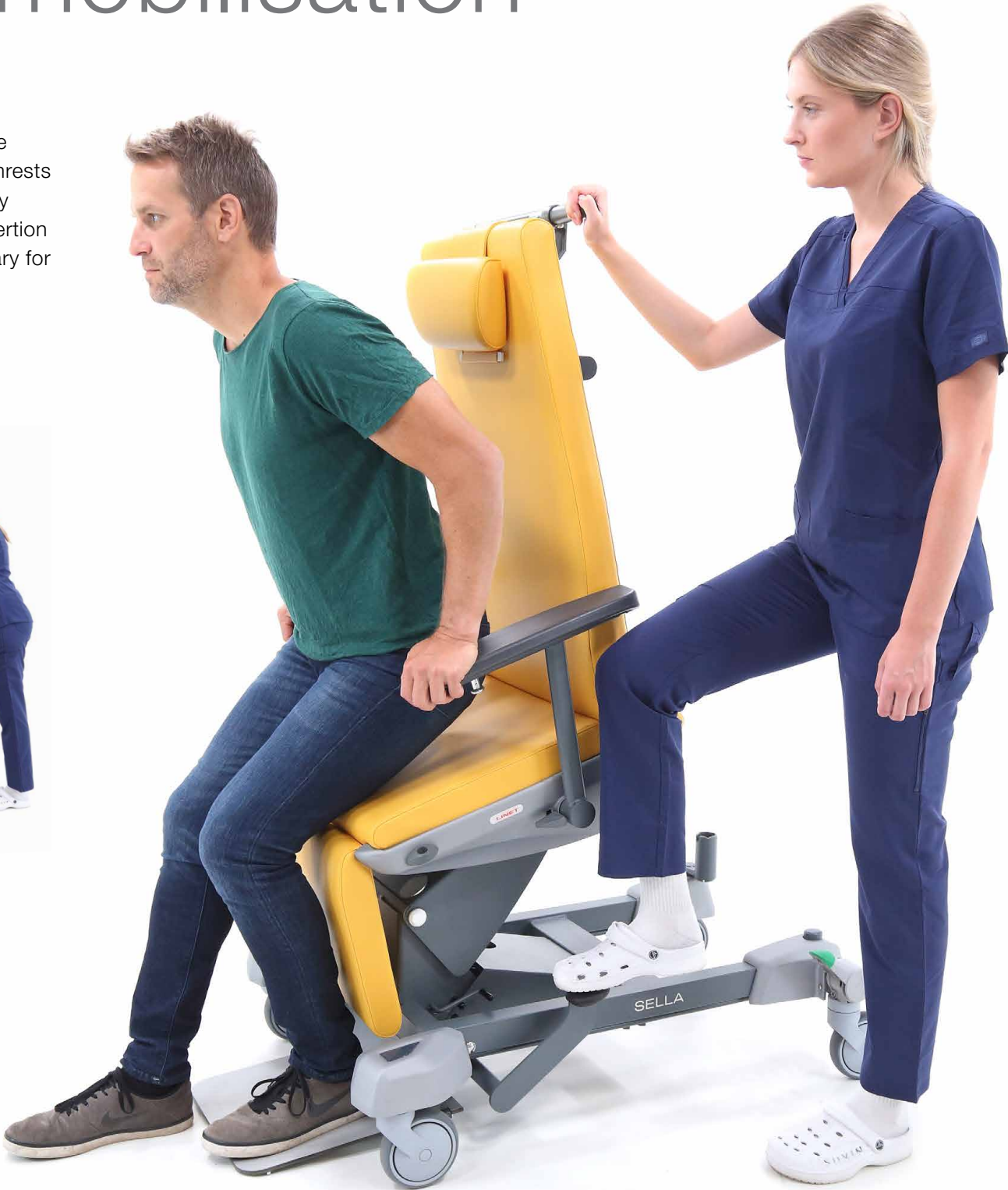

Moreover, the chair brakes can be locked automatically by pressing down on the safety footboard.

Unique variability

High variability of the positions enables patients to adjust the chair according to their individual needs.



Adjustable comfort

Comfortable and easy backrest adjustment enhances overall comfort and supports individual needs.

Positioning for the maximal comfort of the client



Height adjustment



Flat position

For efficient care



Comfortable transfer

Height adjustment makes the work of nurses easier as it enables comfortable transfer of the patient from the chair to the bed.



Easy access

The structure of the Sella chair enables easy access for the staff, for example when ultrasound examination needed.



Simple examination

It is also possible to perform a simple examination or take a blood sample directly from the patient sitting in the chair.

Accessories



Table with safety barrier



Side support (pair)



Stacking pocket



Headrest, big

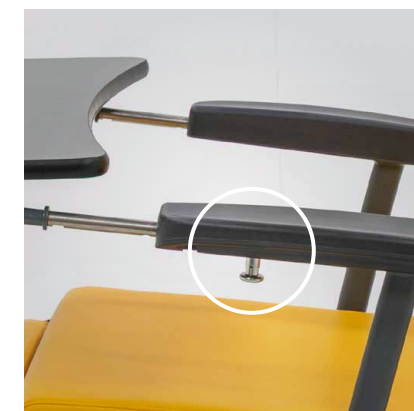


Table safety lock



Urine bag holder



Table holder



Telescopic infusion stand



Holder for infusion stand

Easy cleaning

The construction and used materials of the chair enable easy cleaning without great physical exertion.



Removable seat section for easier maintenance.



Easy and quick cleaning with common cleaning and disinfecting agents.

Colour options

Standard colours



Non standard colours



Technical parameters

	Fixed carriage with individual castor brakes
	Fixed carriage with central castor brakes
	Lift carriage with central castor brakes
	Lift carriage with central castor brakes and help with getting out

Dimensions

Total length in the default position	910 mm
Total length when lying	1,636 mm
Total width	705 mm*
Vertical lift	240 mm
Minimum height***	505 – 510 mm**
Maximum height***	510 – 745 mm**
Back rest adjustment angle	12° to 70°
Footrest adjustment angle	12° to 87°
Arm rest adjustment angle	-25°, 0°, 25°, 90°

Max. loading

Maximum patient weight	150 kg
------------------------	--------

* armrests standard/armrests for blood donation
** by configuration
*** height of the supporting area including upholstery (seat)



LINET

Members of LINET Group

LINET spol. s r.o. | Želevčice 5, 274 01 Slaný | Czech Republic
tel.: +420 312 576 400 | fax: +420 312 522 668 | e-mail: info@linet.com | www.linet.com

